Medicines & Healthcare products Regulatory Agency



Antiserum Reagent 77/568 should be used according to the method described by Wood, JM, Schild, GC, Newman, RW and Seagroatt, VA. Journal of Biological Standardisation, 1977, 5, 237-247.

8. STABILITY

It is the policy of WHO not to assign an expiry date to their international reference materials. They remain valid with the assigned potency and status until withdrawn or amended.

Reference materials are held at NIBSC within assured, temperaturecontrolled storage facilities. Reference Materials should be stored on receipt as indicated on the label. Once reconstituted, diluted or aliquoted, users should determine the stability of the material according to their own method of preparation, storage and use.

NIBSC follows the policy of WHO with respect to its reference materials.

Users who have data supporting any deterioration in the characteristics of any reference preparation are encouraged to contact NIBSC.Reference materials are held at NIBSC within assured, temperature-controlled storage facilities. Reference Materials should be stored on receipt as indicated on the label.

9. REFERENCES

None

10. ACKNOWLEDGEMENTS None

11. FURTHER INFORMATION

Further information can be obtained as follows; This material: enquiries@nibsc.org WHO Biological Standards: http://www.who.int/biologicals/en/ JCTLM Higher order reference materials: http://www.bipm.org/en/committees/jc/jctlm/ Derivation of International Units: http://www.nibsc.org/standardisation/international_standards.aspx Ordering standards from NIBSC: http://www.nibsc.org/products/ordering.aspx NIBSC Terms & Conditions: http://www.nibsc.org/terms_and_conditions.aspx

12. CUSTOMER FEEDBACK

Customers are encouraged to provide feedback on the suitability or use of the material provided or other aspects of our service. Please send any comments to enquiries@nibsc.org

13. CITATION

In all publications, including data sheets, in which this material is referenced, it is important that the preparation's title, its status, the NIBSC code number, and the name and address of NIBSC are cited and cited correctly.

14. MATERIAL SAFETY SHEET

Classification in accordance with Directive 2000/54/EC, Regulation (EC) No 1272/2008: Not applicable or not classified

Physical and Chemical properties			
Physical appearance: White powder	Corrosive:	No	

Influenza Reagent Influenza anti-B/Hong Kong/8/73 HA serum NIBSC code: 77/568 Instructions for use (Version 4.0, Dated 29/03/2008)

1. INTENDED USE

Influenza antiserum reagent 77/568 is prepared in sheep for single radial diffusion assay of B/Hong Kong/8/73 antigens. An appropriate NIBSC reagent should be included in each assay.

2. CAUTION

This preparation is not for administration to humans or animals in the human food chain.

The material is not of human or bovine origin. As with all materials of biological origin, this preparation should be regarded as potentially hazardous to health. It should be used and discarded according to your own laboratory's safety procedures. Such safety procedures should include the wearing of protective gloves and avoiding the generation of aerosols. Care should be exercised in opening ampoules or vials, to avoid cuts.

3. UNITAGE

No unitage is assigned to this material.

4. CONTENTS

Country of origin of biological material: United Kingdom. 77/568 antiserum reagent is freeze dried sheep serum.


5. STORAGE

+2-8⁰C

Please note: because of the inherent stability of lyophilized material, NIBSC may ship these materials at ambient temperature.

6. DIRECTIONS FOR OPENING

Tap the ampoule gently to collect the material at the bottom (labelled) end. Ensure ampoule is scored all round at the narrow part of the neck, with a diamond or tungsten carbide tipped glass knife file or other suitable implement before attempting to open. Place the ampoule in the ampoule opener, positioning the score at position 'A'; shown in the diagram below. Surround the ampoule with cloth or layers of tissue paper. Grip the ampoule and holder in the hand and squeeze at point 'B'. The ampoule will snap open. Take care to avoid cuts and projectile glass fragments that enter eyes. Take care that no material is lost from the ampoule and that no glass falls into the ampoule.

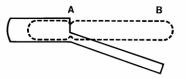

Side view of ampoule opening device containing an ampoule positioned ready to open. 'A' is the score mark and 'B' the point of applied pressure.

7. USE OF MATERIAL

No attempt should be made to weigh out any portion of the material

For the assay of antigens containing 50-100 μ g of HA activity in 1ml, approximately 25 μ l of the undiluted Reagent should be added to 1ml agarose. It may be necessary to change the antiserum concentrations according to local laboratory conditions.

National Institute for Biological Standards and Control,

Potters Bar, Hertfordshire, EN6 3QG. T +44 (0)1707 641000, nibsc.org WHO International Laboratory for Biological Standards,

UK Official Medicines Control Laboratory

Medicines & Healthcare products Regulatory Agency



Stable: Yes	Oxidising:	No		
	luuite et.	Nia		
Hygroscopic:	Irritant:	No		
No				
Flammable:	Handling:	See caution, Section 2		
No				
Other (specify): Contains Sheep Serum and Sodium Azide (0.05% w/v)Please complete				
Toxicological properties				
Effects of inhalation: Avoid inhalation				
Effects of ingestion: Avoid ingestion				
Effects of skin absorption: Avoid contact with skin				
Suggested First Aid				
Inhalation: S	eek medical advice			
Ingestion: Seek medical advice				
Contact with eyes: W medical advice	ash with copious	s amounts of water. Seek		
Contact with skin: W	ash thoroughly v	vith water.		
Action on Spillage and Method of Disposal				
material wetted with a appropriate disinfecta	n appropriate dis nt followed by wa	e taken up with absorbent infectant. Rinse area with an iter. age should be treated as		

Absorbent materials used to treat spillage should be treated as biological waste.

15. LIABILITY AND LOSS

In the event that this document is translated into another language, the English language version shall prevail in the event of any inconsistencies between the documents.

Unless expressly stated otherwise by NIBSC, NIBSC's Standard Terms and Conditions for the Supply of Materials (available at http://www.nibsc.org/About_Us/Terms_and_Conditions.aspx or upon request by the Recipient) ("Conditions") apply to the exclusion of all other terms and are hereby incorporated into this document by reference. The Recipient's attention is drawn in particular to the provisions of clause 11 of the Conditions.

16. INFORMATION FOR CUSTOMS USE ONLY Country of origin for customs purposes*: United Kingdom * Defined as the source update the goade have been produced

* Defined as the country where the goods have been produced and/or sufficiently processed to be classed as originating from the country of supply, for example a change of state such as freeze-drying. Net weight: 1g Toxicity Statement: Non toxic

Veterinary certificate or other statement if applicable. Attached: No

National Institute for Biological Standards and Control, Potters Bar, Hertfordshire, EN6 3QG. T +44 (0)1707 641000, nibsc.org WHO International Laboratory for Biological Standards, UK Official Medicines Control Laboratory