



**Non WHO Reference Material  
RENIN (rDNA, Human Type)  
NIBSC code: 89/544  
Instructions for use  
(Version 5.0, Dated 30/04/2013)**

**This material is not for in vitro diagnostic use.**

**1. INTENDED USE**

Calibrant for renin assays

**2. CAUTION**

**This preparation is not for administration to humans or animals in the human food chain.**

The preparation contains material of human origin, and either the final product or the source materials, from which it is derived, have been tested and found negative for HBsAg, anti-HIV and HCV RNA. As with all materials of biological origin, this preparation should be regarded as potentially hazardous to health. It should be used and discarded according to your own laboratory's safety procedures. Such safety procedures should include the wearing of protective gloves and avoiding the generation of aerosols. Care should be exercised in opening ampoules or vials, to avoid cuts.

**3. UNITAGE**

To be calibrated in terms of International Units (IU) of renin using the IS for renin (preparation 68/356).

**4. CONTENTS**

Country of origin of biological material: Renin: USA, HSA: UK.

Each ampoule contains the residue after freeze-drying of 0.5ml of a solution that contained:-

1.53 µg	Renin (rDNA human type)
3.0 mg	Trehalose
1.0 mg	human serum albumin
0.05 mol/l	sodium phosphate buffer, pH 6.5

The material has not been sterilized and the ampoules contain no bacteriostat.


**5. STORAGE**

Store at -20°C or below

**Please note: because of the inherent stability of lyophilized material, NIBSC may ship these materials at ambient temperature.**

**6. DIRECTIONS FOR OPENING**

Tap the ampoule gently to collect the material at the bottom (labelled) end. Ensure ampoule is scored all round at the narrow part of the neck, with a diamond or tungsten carbide tipped glass knife file or other suitable implement before attempting to open. Place the ampoule in the ampoule opener, positioning the score at position 'A'; shown in the diagram below. Surround the ampoule with cloth or layers of tissue paper. Grip the ampoule and holder in the hand and squeeze at point 'B'. The ampoule will snap open. Take care to avoid cuts and projectile glass fragments that enter eyes. Take care that no material is lost from the ampoule and that no glass falls into the ampoule.

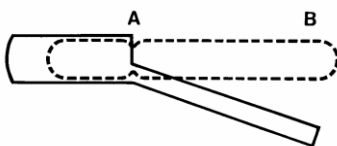

Side view of ampoule opening device containing an ampoule positioned ready to open. 'A' is the score mark and 'B' the point of applied pressure.

**7. USE OF MATERIAL**

No attempt should be made to weigh out any portion of the freeze-dried material prior to reconstitution.

For all practical purposes each ampoule contains the same amount of the same materials. Dissolve the total contents in a known amount of suitable buffer solution with carrier protein (free of peptidase), where extensive dilution is required, to minimise loss by surface adsorption.

For economy of use it is recommended that the solution be subdivided into several small containers, which are frozen rapidly to below -40°C and then stored below -40°C in the dark. Repeated freezing and thawing should be avoided.

**8. STABILITY**

NIBSC follows the policy of WHO with respect to its reference materials. It is the policy of WHO not to assign an expiry date to their international reference materials. They remain valid with the assigned potency and status until withdrawn or amended. Reference materials are held at NIBSC within assured, temperature-controlled storage facilities. Reference Materials should be stored on receipt as indicated on the label. In addition, once reconstituted, diluted or aliquoted, users should determine the stability of the material according to their own method of preparation, storage and use. Users who have data supporting any deterioration in the characteristics of any reference preparation are encouraged to contact NIBSC.

**9. REFERENCES**

Not applicable

**10. ACKNOWLEDGEMENTS**

The renin contained in this ampoule was selected from preparations generously contributed to WHO by the Upjohn Company and CALBIO.

**11. FURTHER INFORMATION**

Further information can be obtained as follows:

This material: [enquiries@nibsc.org](mailto:enquiries@nibsc.org)

WHO Biological Standards:

<http://www.who.int/biologicals/en/>

JCTLM Higher order reference materials:

<http://www.bipm.org/en/committees/jc/jctlm/>

Derivation of International Units:

[http://www.nibsc.org/standardisation/international\\_standards.aspx](http://www.nibsc.org/standardisation/international_standards.aspx)

Ordering standards from NIBSC:

<http://www.nibsc.org/products/ordering.aspx>

NIBSC Terms & Conditions:

[http://www.nibsc.org/terms\\_and\\_conditions.aspx](http://www.nibsc.org/terms_and_conditions.aspx)

**12. CUSTOMER FEEDBACK**

Customers are encouraged to provide feedback on the suitability or use of the material provided or other aspects of our service. Please send any comments to [enquiries@nibsc.org](mailto:enquiries@nibsc.org)

**13. CITATION**

In all publications, including data sheets, in which this material is referenced, it is important that the preparation's title, its status, the NIBSC code number, and the name and address of NIBSC are cited and cited correctly.

**14. MATERIAL SAFETY SHEET**

Classification in accordance with Directive 2000/54/EC, Regulation (EC) No 1272/2008: Not applicable or not classified

<b>Physical and Chemical properties</b>	
Physical appearance: Freeze-dried powder	Corrosive: No
Stable: Yes	Oxidising: No
Hygroscopic: Yes	Irritant: No
Flammable: No	Handling: See caution, Section 2
Other (specify): Contains material of human origin	
<b>Toxicological properties</b>	
Effects of inhalation:	Not established, avoid inhalation
Effects of ingestion:	Not established, avoid ingestion
Effects of skin absorption:	Not established, avoid contact with skin
<b>Suggested First Aid</b>	
Inhalation:	Seek medical advice
Ingestion:	Seek medical advice
Contact with eyes:	Wash with copious amounts of water. Seek medical advice
Contact with skin:	Wash thoroughly with water.
<b>Action on Spillage and Method of Disposal</b>	
Spillage of contents should be taken up with absorbent material wetted with an appropriate virucidal agent. Rinse area with an appropriate virucidal agent followed by water. Absorbent materials used to treat spillage should be treated as biological waste.	

#### 15. LIABILITY AND LOSS

In the event that this document is translated into another language, the English language version shall prevail in the event of any inconsistencies between the documents.

Unless expressly stated otherwise by NIBSC, NIBSC's Standard Terms and Conditions for the Supply of Materials (available at [http://www.nibsc.org/About\\_Us/Terms\\_and\\_Conditions.aspx](http://www.nibsc.org/About_Us/Terms_and_Conditions.aspx) or upon request by the Recipient) ("Conditions") apply to the exclusion of all other terms and are hereby incorporated into this document by reference. The Recipient's attention is drawn in particular to the provisions of clause 11 of the Conditions.

#### 16. INFORMATION FOR CUSTOMS USE ONLY

<b>Country of origin for customs purposes*:</b> United Kingdom * Defined as the country where the goods have been produced and/or sufficiently processed to be classed as originating from the country of supply, for example a change of state such as freeze-drying.
<b>Net weight:</b> 4mg
<b>Toxicity Statement:</b> Toxicity not assessed
<b>Veterinary certificate or other statement</b> if applicable. <b>Attached:</b> No